
DATA SHEET – Mattress MTH-01



The MTH01 heat-sealed mattress is designed for secure healthcare and psychiatric environments. Manufactured in France, it combines high-resilience polyurethane foam with a 100% waterproof PVC cover to ensure comfort, hygiene, and durability for intensive use.

Thanks to its seamless heat-sealed construction, the MTH01 mattress provides optimal waterproofing and excellent tear resistance. Compatible with hospital disinfection protocols, it also offers outstanding durability, easy maintenance, and enhanced protection against fluid ingress and contamination.

▪ **Dimensions**	L 200 cm x W 90 cm L 200 cm x W 120 cm
▪ **Thickness**	20 cm
▪ **Weight**	4 kg
▪ **Options**	Blue

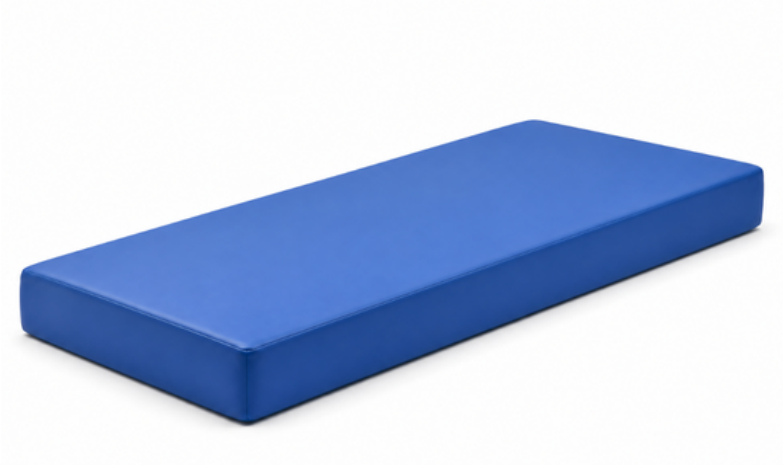

Click here to learn more

Standards and compliance

- **REACH (EC No. 1907/2006)** – Regulation on chemicals and chemical substances
- **ISO 846:2019** – Resistance to fungal and bacterial growth
- **NF EN 597** – Assessment of the ignitability of mattresses
- Manufactured under **MASE certification**
- **ISO 9001** certification process underway

Cleaning Protocol

- Washable and waterproof upholstery, suitable for hospital environments
- Compatible with hospital cleaning and disinfection protocols
- Pre-cleaning: remove visible dirt and contaminants before disinfection
- Routine cleaning with a neutral detergent, using a soft, non-abrasive cloth or sponge
- Thoroughly clean the entire surface, particularly seams, welded joints, and connecting areas
- Disinfect using a product compatible with the mattress upholstery, following the manufacturer's recommended concentrations and contact times
- Do not use abrasive products, solvents, or chemicals that may damage the upholstery
- Rinse if necessary, according to the cleaning product manufacturer's instructions
- Air dry or wipe with a clean, lint-free cloth
- Allow the mattress to dry completely before reuse
- Regular cleaning recommended in hospital environments
- Clean and disinfect between patients according to the facility's hygiene protocols
- Regularly inspect the upholstery for signs of damage, including cracks, tears, punctures, or damaged seams
- If damage compromises the waterproof integrity of the upholstery, remove the mattress from service